





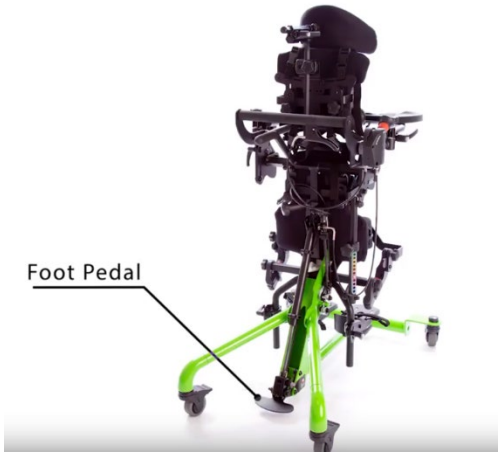



# EasyStand Bantam Extra Small and Small Mobile Moderate Support Package Feature Descriptions

FEATURE	IMAGE	DESCRIPTION
<p><b>Tall Mobile Frame</b></p>		<p>This feature allows self-propulsion. Push rims adjust forward and aft. Mobile footprint 24"x38" (61 x 97cm). Push rim height 35" (89cm) from floor. 7" Front Wheels are standard forward and backward rolling wheels.</p>
<p><b>Multi-Adjustable Foot Plates</b></p>		<p>Foot plates are adjustable forward/aft +/- .75" (2cm), plantar/dorsi +/- 20°, toe-in/toe-out, swivel 30° left/right, in addition to height.</p>
<p><b>Secure Foot Straps</b></p>		<p>Velcro® and D-ring adjustment holds feet in alignment with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro overlap. Must choose size in drop box.</p>

FEATURE	IMAGE	DESCRIPTION
<b>Swing-Away Knee Pads</b>		<p>The push of a button swings knee pad away providing more room for transfers. Multi-adjustable knee pads are adjustable in swivel, height +/- .5" (+/-1.3cm), width between knee pads (measured from the center of one knee pad to the center of the other knee pad) 5.5"-10" (14-25cm) in addition to a depth adjustment range of 2.5"-8" (6-20cm). Must choose inside width of knee pads in drop box.</p>
<b>Planar Seat and Back</b>		<p>Flat Seat and Back provide additional support. Must choose size in drop box. These measurements represent the distance from the seat to the top of the back and are given in a range because the seat back height is adjustable.</p>
<b>Easy Adjust Seat Depth</b>		<p>Roller bearing seat depth allows for smoother seat depth adjustment.</p>
<b>Velcro® Positioning Belt</b>		<p>2" (5cm) wide positioning belt has a D-ring Velcro closure, provides hip stability for the user. Must choose size in drop box.</p>

FEATURE	IMAGE	DESCRIPTION
<b>Hip Supports</b>		<p>Supports the hips in symmetrical alignment from sitting to standing. Adjustable width from 6"-9" (15-23cm) on extra small models and 8"-12" (20-30cm) on small models. 3 settings with a 2" (5cm) height range. Pad size is 3.5"x3.5" (9x9cm).</p>
<b>Chest Strap</b>		<p>Velcro® and D-ring Adjustment Chest strap provides moderate support to help maintain the user's trunk in midline.</p>
<b>X-Style Chest Vest with High Mount Chest Vest Bracket</b>		<p>Neoprene® with camlock adjustment, the X-Style Chest Vest provides anterior trunk support for the child through four points of control. Support is provided to help maintain the user's trunk in midline and shoulders back. Must choose size in drop box. The High Mount Chest Vest Bracket moves mounting points for chest vest higher to secure shoulders in position.</p>
<b>Mobile Tray</b>		<p>Tray provides anterior support while using the mobile option. Tray is 4"x12" (10x30cm) and has height and depth adjustment.</p>

FEATURE	IMAGE	DESCRIPTION
<p><b>Gas Spring Lift with Foot Pedal</b></p>		<p>The Gas Spring Lift Mechanism is used to adjust the angle and the position of the stander and is controlled by the Foot Pedal.</p>
<p><b>Swivel Castors</b></p>		<p>Locking castors allow for easier steering and movement of the stander.</p>